

Requesting Add-On Tests

Easily request add-on tests for existing orders directly from your electronic health record (EHR). Enabled users can access add-on test requests from the original result. Add-on tests cannot be requested for open orders.

Note: To enable add-on test requests in Diagnostic Assistant, users must have a Labcorp Link™ account.

1. Launch the Add-On Testing tab

A Click the **Request Add-On** icon.

Note: The first time a user accesses a Labcorp Link-enabled feature, such as add-on tests, they'll be asked to sign in to Labcorp Link and allow access. Their credentials will be used to enable these features as part of a **one-time** authorization.

The screenshot shows the 'Laboratory Results' interface with a table of test results. A callout box highlights the 'Request Add-On' icon in the 'Actions' column. Another callout box shows the Labcorp sign-in process, including a 'Sign In' button and an 'Allow Access' button.

Date Collected	Test Name	Status	Source	Actions
03/22/2024 10:45 AM	Hemoglobin A1c	Final	Labcorp	[Icon] +
03/18/2024 12:00 AM	Specimen Status	Registered	Labcorp	[Icon] +
03/18/2024 12:00 AM	Ca+PTH Intact	Final	Labcorp	[Icon] +
03/18/2024 12:00 AM	Hemoglobin A1c	Final	Labcorp	[Icon] +
03/11/2024 12:00 AM	Celiac Ab tTG TIGa w/RtX	Final	Labcorp	[Icon] +
03/04/2024 10:46 AM	Thyroid Panel With TSH	Final	Labcorp	[Icon] +
03/04/2024 10:15 AM	Hemoglobin A1c	Final	Labcorp	[Icon] +
02/29/2024 12:00 AM	Celiac Ab tTG TIGa w/RtX	Final	Labcorp	[Icon] +
02/22/2024 05:45 AM	Celiac Ab tTG TIGa w/RtX	Final	Labcorp	[Icon] +
02/16/2024 10:00 AM	Celiac Ab tTG TIGa w/RtX	Final	Labcorp	[Icon] +

2. Review the original order information

A Click **View Details** to display existing add-on requests for the order. Duplicate test requests on the same order are not permitted.

The screenshot shows the 'Add-On Testing' tab in the Laboratory Results interface. A callout box highlights the 'View Details' link in the 'Add On Requests (1):' section. Another callout box shows the 'Add-On Request Details' dialog, which includes a table of external requests.

Order Date	Ordered Tests	Request Status	Accession
Not Available	005058 - Hematocrit 977900 - Written Authorization	Some results have been entered	0624802001

Requesting Add-On Tests

3. Enter add-on test information

- A** Type a test number or name in the **Add Test** field.
- B** To select from a picklist, click the **All Tests** drop-down and choose the desired picklist. Then, select from tests included in the picklist.

Note: A result cannot have more than two add-on tests. Specimen requirements for the added tests must match the specimens received for the original order. Other restrictions apply.

Laboratory Results Maria Cruz
03/09/1994

PanelView ClinicalView Add-On Testing x

Add-On Tests

Add Test * **All Tests** v

- 001032 - Glucose
- 090399 - Glucose Tolerance (7 Sp Blood)
- 090407 - Glucose Tolerance (8 Sp Blood)
- 090423 - Glucose/Insulin Resp (7 Spec)
- 211508 - Glucose Tolerance Prof (7 Sp)
- 001818 - Glucose, Plasma
- 211516 - Glucose Tolerance Prof (8 Sp)
- 239822 - Glucose (2 Spec) Tolerance, S

Selected Tests

001032 - Glucose x

Reminder: When ordering tests for which Medicare or Medicaid reimbursement is sought, the provider should only order those tests which the provider believes are medically necessary for each patient.

Add-On Diagnoses

Add Diagnosis **All Diagnoses** v

Selected Diagnoses

N20.0 - CALCULUS OF KIDNEY (Original Diagnosis)

The ordering physician must always determine, for the specific date of service, the

4. Enter diagnosis information (optional)

- A** Type a diagnosis code or name in the **Add Diagnosis** field. If no selection is made, the diagnosis code from the original order will be used.
- B** To select from a picklist, click the **All Diagnoses** drop-down and choose the desired picklist. Then, select from tests included in the picklist.

Laboratory Results Maria Cruz
03/09/1994

PanelView ClinicalView Add-On Testing x

Add-On Diagnoses

Add Diagnosis **All Diagnoses** v

- G93.31 - POSTVIRAL FATIGUE SYNDROME
- G93.32 - MYALGIC ENCPHLOMYEL/CHR FATIGUE SYN
- G93.39 - OTH POST INFECT&RELATED FATIGUE SYN
- M48.40XA - FATIGUE FX VERT SITE UNS INIT ENC
- M48.40XD - FATIGUE FX VERT UNS SUB ROUTIN HEAL
- M48.40XG - FATIGUE FX VERT UNS SUB DELAY HEAL
- M48.40XS - FATIGUE FX VERT SITE UNS SEQUELA FX
- M48.41XA - FATIGUE FX VERT OAA REGION INIT ENC

Selected Diagnoses

N20.0 - CALCULUS OF KIDNEY (Original Diagnosis)

G93.31 - POSTVIRAL FATIGUE SYNDROME x

The ordering physician must always determine, for the specific date of service, the appropriate diagnosis codes based on the patient's signs and symptoms.

Point Of Contact

Who should we contact if we have questions about this request?

Requesting Add-On Tests

5. Confirm requester information and submit the request

When an add-on test request is submitted, a Labcorp customer service agent reviews it to determine the specimen type required, specimen age and the remaining volume. To edit or cancel an add-on test request, contact Labcorp at 855-LABCORP (855-522-2677).

Results for add-on tests appear in the EHR and in Diagnostic Assistant.*

*Unsolicited results: For orders originally placed in an electronic health record (EHR), add-on requests may be considered "unsolicited results" and routed to an error queue for special handling as results file back in. Check your EHR specifications to see how this may work for you.

Laboratory Results Maria Cruz
03/09/1994

Panel View Clinical View **Add-On Testing**

sought, the provider should only order those tests which the provider believes are medically necessary for each patient.

Add-On Diagnoses

Add Diagnosis

Enter Diagnoses Code or Name All Diagnoses

Selected Diagnoses

- N20.0 - CALCULUS OF KIDNEY (Original Diagnosis)
- G93.31 - POSTVIRAL FATIGUE SYNDROME

The ordering physician must always determine, for the specific date of service, the appropriate diagnosis codes based on the patient's signs and symptoms.

Point Of Contact

Who should we contact if we have questions about this request?

Name	Phone Number
<input type="text" value="User Name"/>	<input type="text"/>


Note: Results for added tests may not file into EHR systems that do not accept unsolicited results.

Submit Request Cancel

Requesting Add-On Tests

Specimen Requirements: When adding tests, specimen requirements for the added tests must match the specimens received for the original order.

- A** Click the **How to Resolve** link for a mismatched test to learn more on how to fix the issue.
- B** Click the test name to review specimen requirements for a selected test.

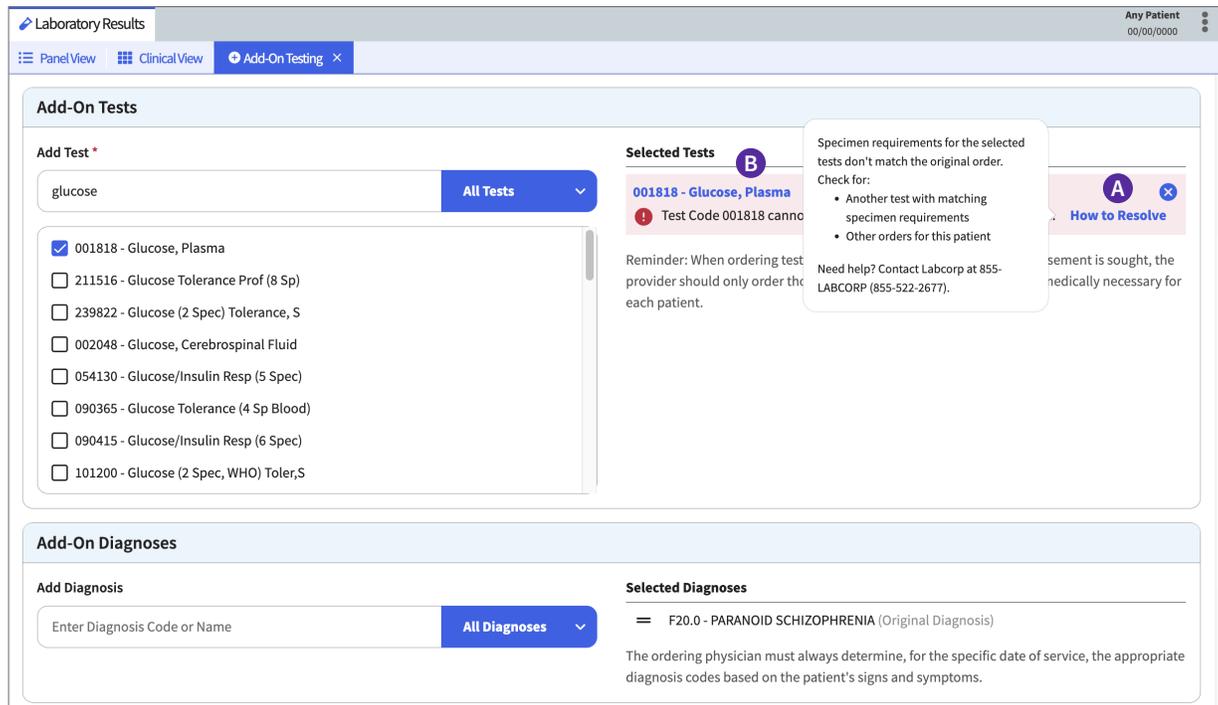

The screenshot displays the 'Add-On Testing' interface in the Labcorp Laboratory Results system. The 'Add Test' field is set to 'glucose'. A list of available tests is shown, with '001818 - Glucose, Plasma' selected. The 'Selected Tests' section shows '001818 - Glucose, Plasma' with a red error icon and a 'How to Resolve' link. A tooltip explains that specimen requirements for the selected tests don't match the original order and provides a list of checks: 'Another test with matching specimen requirements' and 'Other orders for this patient'. A reminder states: 'When ordering test, provider should only order the test for each patient. Need help? Contact Labcorp at 855-LABCORP (855-522-2677).'

Additional Requirements: Certain criteria must be met for the Request Add-On icon to appear.

The order must:

- Be the original Labcorp order
- Be resulted, with an accession number ending in 0
- Have a collection date within the last 8 days

The user must have:

- An Administrator, Clinical or Provider role within Labcorp Link (i.e., not a Guest user)
- Results access
- Access to the ordering account of the original order



For more information,
contact your local Labcorp representative.